



# PATELLA FRACTURE REPAIR & SOFT-TISSUE MANAGEMENT, RIGHT OUT OF THE BOX.

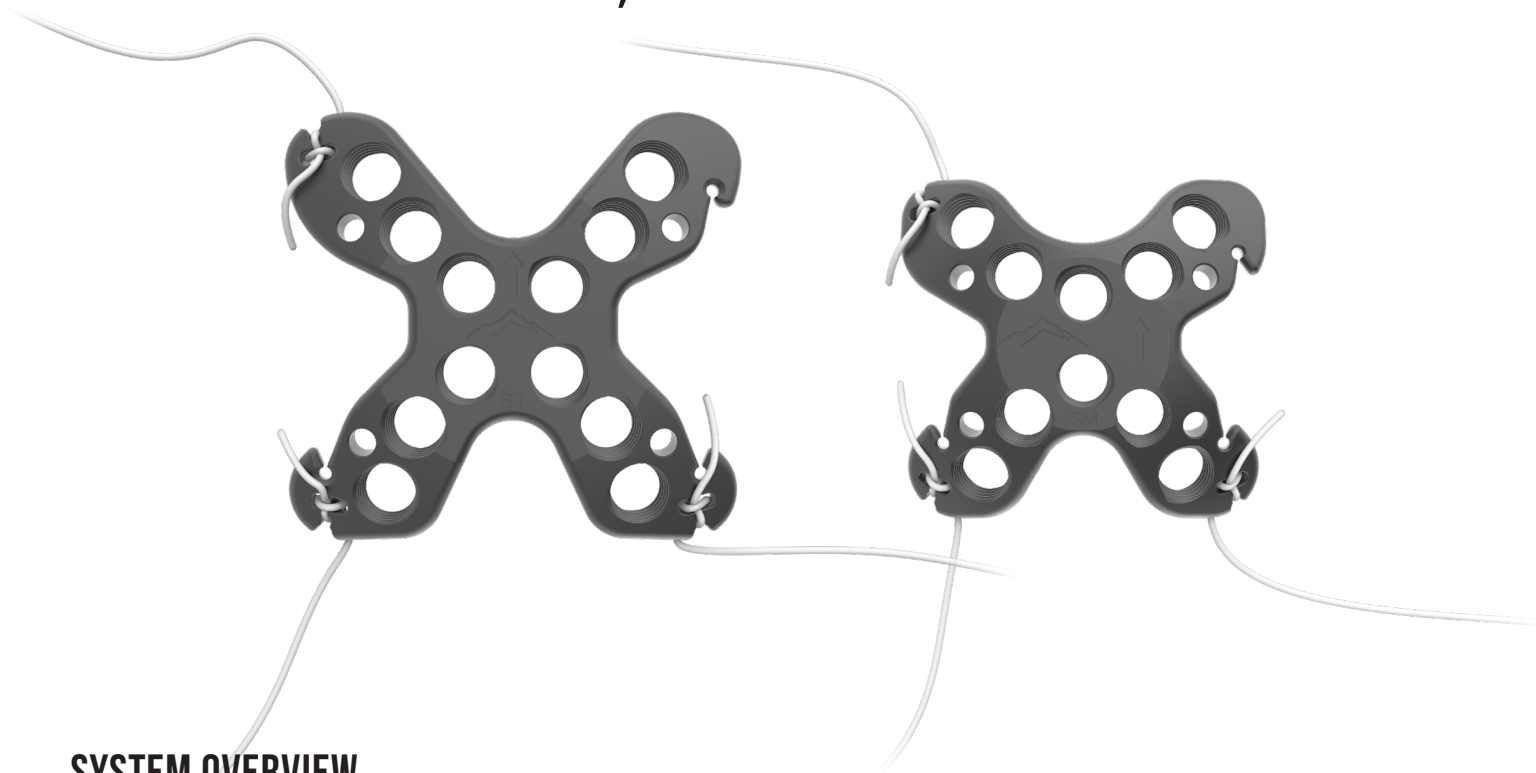

## SYSTEM OVERVIEW

The Summit Patella Plating System is a sterile, single-use procedure kit indicated to treat patella fractures in adults. The kit is available with either a Standard or Small sized plate. Each kit comes ready to use with a pre-loaded #2 suture for the quadriceps tendon and two pre-loaded #2 sutures for the patellar ligament. The sutures in the inferior pole are designed to allow for an optional Krackow stitch technique.

Once the sutures have been placed, the surgeon can cinch them tight by simply pulling on the tail of the suture toward the center of the plate.

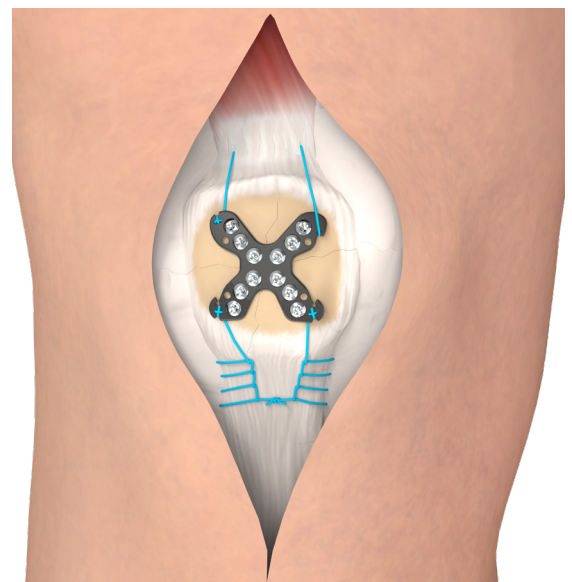

*Note: Blue sutures are for visual reference only. Sutures included in procedure kits are standard white.*

# SUMMIT PATELLA PLATING SYSTEM



## SYSTEM COMPONENTS

- 1 LOCKING DRILL GUIDE**  
Locking drill guide with set angle
- 2 T8 DRIVER SHAFT**  
AO Quick Connect T8 driver shaft
- 3 1.8 MM DRILL**  
1.8 mm AO Quick Connect solid drill bit
- 4 OLIVE WIRE**  
Provides provisional fixation
- 5 1.6 MM K-WIRE**  
Provides provisional fixation and fracture reduction
- 6 T8 DRIVER**  
Driver with handle
- 7 VA CONE DRILL GUIDE**  
Allows for drilling up to 15° off axis in any direction
- 8 LOCKING SCREWS**  
2.4 mm locking screws available in 12 mm\*, 14 mm, 16 mm, and 18 mm\* lengths
- 9 PATELLA PLATE**  
Patella-specific plate design with integrated sutures, available in small and standard sizes
- 10 SUTURES**  
A single #2 suture for the quadriceps tendon and two #2 sutures for the patella ligament

\* 12 mm - available with Small Plate Kit only

\* 18 mm - available with Standard Plate Kit only



## PRODUCT SPECIFICATIONS

### PLATE DIMENSIONS

#### Standard Plate:

Length: 30.5 mm  
Width: 29.1 mm  
Thickness: 1.25 mm

#### Small Plate:

Length: 26 mm  
Width: 25 mm  
Thickness: 1.4 mm

### PART NUMBERS:

Standard Plate Kit: **1316-3000**

Small Plate Kit: **1316-3001**

## CONTACT

U.S. Toll Free:  
800.818.4726

[endeavorortho.com](http://endeavorortho.com)